


Net gas type safety cabinet

PRODUCT INTRODUCTION

SF-Z filtration series high performance ductless filtered safety storage cabinets are produced strictly in accordance with OSHA29 CFR 1910.106, NFPA CODE30 and Q/320205 XBB001-2014 standard, and pass through a third party authority the European CE certification and ISO9001 certification. The quality of product can meet and exceed the requirement of international market for safety storage product.

EXCELLENT

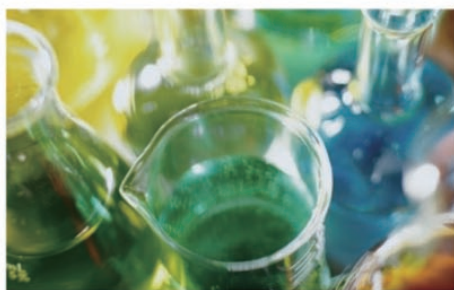


- ◆ No duct, easy to install, no exhaust gas, modern and environment protection
- ◆ Optional filtration module system according to a variety of experiment needs
- ◆ Advanced module filtration technology, full absorption of toxic vapor, particle & dust, etc.
- ◆ No consumption of air conditioning energy, high efficiently saving energy
- ◆ Moving conveniently, nearby storage, easy to access, to improve the work efficiency

WORKING PRINCIPLE

③ Choose a proper high efficient filter according to the chemicals

④ The hazardous air goes through molecular filter become fresh air.

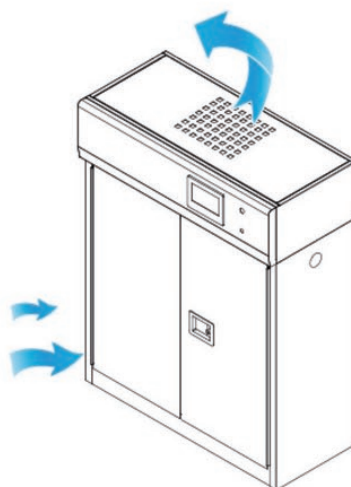

① The turbine motor draws the air from the outside of cabinet.

② The air brings the hazardous vapour into the filtration system.

NET GAS TYPE SAFETY CABINET-product feature

Activated carbon filter system with
US imported fan

Touch screen computer control system real
time monitor the humidity temperature and
concentration in the cabinet




Reflective warning label

Galvanized removable shelf

Double wall 1.2mm steel plate
38mm air insulating space



Vent with anti-prevention

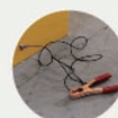

accommplate with OSHA
standard static grounding
connector

50mm leakage proof sump



NET GAS TYPE SAFETY CABINET-product feature



3-Point locking system:

Double locks +3 latch, easy to open 180 degrees cupboard door, equipped with double key; optional padlock, realize double lock management.



Temperature, humidity VOC alarm device:

- ① A temperature and humidity controller, to monitor the relative temperature and humidity in the cabinet, digital display settings and measured values, it will alarm immediately once the temperature and humidity exceed the setting measured values.
- ② Silent turbine fan draws chemical vapors from bottles into filter system with high absorption capacity. No static, no spark and super silence motor
- ③ Various type filters to realize 99.99% high efficiently filtering.



Net gas type safety cabinet

Model	Specification	Configure	Optional parts
SF-Z030	Dimension(mm): 1400×1090×460mm	Filters(pcs): 4	Shelf Wheel Filter saturation alarm Filters
		Fans(pcs): 1	
	Door: Double door	Display screen: YES	
	Lock: 3-points lock	Power cord(pcs): 1	
	Shelf: 1 galvanized shelf	Control system(set): 1	
	Colour: Yellow	Alarm system: YES	
	Noise: 40dBA		
	Rated Voltage: 220V/50Hz		
	Electric Powder: 40W		
SF-Z045	Rated current: 2A		Shelf Wheel Filter saturation alarm Filters
	Dimension(mm): 1930×1090×460mm	Filters(pcs): 4	
		Fans(pcs): 1	
	Door: Double door	Display screen: YES	
	Lock: 3-points lock	Power cord(pcs): 1	
	Shelf: 1 galvanized shelf	Control system(set): 1	
	Colour: Yellow	Alarm system: YES	
	Noise: 40dBA		
	Rated Voltage: 220V/50Hz		
	Electric Powder: 40W		
	Rated current: 2A		